

# DRUGKEEPER

## Storage Cabinet



For the storage of illegal and hazardous chemicals!

Exceeds OSHA, ANSI and all relevant international standards

### APPLICATIONS

- Protection from 'off-gassing' of typical chemicals found at Meth labs including Methamphetamine
- Protection from odors and particulates while drying marijuana including hazardous mold (aflatoxin)



DRUGKEEPER 34S

### FEATURES

- HEPA and Carbon filters
- Door locks for security
- Adjustable and removable shelving
- Non-transparent tinted door option to prevent viewing inside the cabinet
- Floor standing, bench top, or wall mountable units available



DRUGKEEPER 64T

Shown with clear doors. Optional non-transparent doors available.

Toll Free/ 800.306.0656 F/ 800.306.0677

[WWW.AIR-SCIENCE.COM](http://WWW.AIR-SCIENCE.COM)



# AIR SCIENCE USA



Made in the USA

## DRUGKEEPER STORAGE CABINET

The *DRUGKEEPER* is a containment cabinet for the drying and storing of drugs. The airflow through the cabinet contains any fumes, odors, or particulates given off and ensures a safe working environment for personnel.

### METH LAB CHEMICALS

As the methamphetamine lab problem continues to grow, agencies need to assess the health and environmental risks of handling and storing these types of chemicals during seizure until they can be properly processed and disposed.


Exposure to meth residues may cause symptoms similar to those experienced by meth users. In addition, exposure to commonly used ingredients such as volatile compounds (VOC's) can cause nose irritation, headaches, dizziness, nausea, vomiting, and breathing difficulties. Some of these compounds are known to cause cancer. Acids that are used can cause a burning sensation and severe eye damage.



Interior Shelving

### STORAGE OF MARIJUANA

Exposure to toxic mold on marijuana plants, called aflatoxins, can be extremely dangerous. The *DRUGKEEPER* is equipped with HEPA filters that can contain these molds. In addition the airflow through out the cabinet will help quickly air dry 'wet' plants while the carbon filters within the cabinet will help eliminate the strong odors that are present.

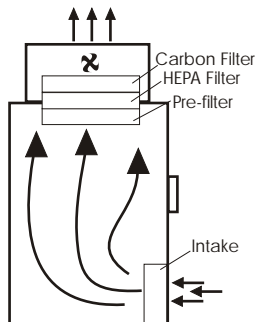

Air Flow Pattern

### STORAGE OF COCAINE & OTHER SUBSTANCES

Most illegal drugs can be safely stored within the *DRUGKEEPER* to prevent operator exposure during processing and disposal.

### FILTERS

Each *DRUGKEEPER* is equipped with a pre-filter, high capacity carbon filters to contain chemical vapors and odors, and HEPA filters to contain particulates and mold.

### FEATURES

- All chemical resistant polypropylene construction
- Free standing or bench top designs
- Lockable tinted front door with warning stickers
- Bottom spill tray to contain potential spills
- Fully adjustable shelving
- Adjustable leveling feet or locking wheels
- Locking door to restrict unauthorized entry
- Low airflow alarm that alerts operator when airflow drops to an unacceptable level or in case of fan failure.



Lockable Wheels

For advice on filter selection call Air Science USA 800.306.0656



610 Center Road • PO Box 62296 • Fort Myers, Florida 33906 • T/239.489.0024 • F/800.306.0677 • Toll Free/800.306.0656 • www.air-science.com



### SPECIFICATIONS ON ALL UNITS

Airflow (CFM)	175
Filters	HEPA and activated carbon
Prefilter (Electrostatic)	Filter efficiencies superior to 99.6%
Supply (single phase)	110V 60Hz single phase (Other voltages available)
Fan	Centrifugal brushless & sparkless
Alarms	Low airflow alarm
Noise (dBA)	<49
Cabinet	White polypropylene
Door Mechanism	Clear or Solid with key lock

### DIMENSIONS

Drugkeeper 34S	(")	34w x 20d x 29h
Drugkeeper 34T	(")	34w x 24d x 79h
Drugkeeper 64T	(")	64w x 24d x 79h

### NUMBER OF SHELVES

Drugkeeper 34S	4
Drugkeeper 34T	5
Drugkeeper 64T	10

Custom Sizes Available

Distributed by: